



PORTABLE BARIATRIC FLOOR SCALE 6800

24-INCH-SQUARE LOW-PROFILE
PLATFORM WITH BUILT-IN
CARRYING HANDLE



Battery or AC Power • BMI Calculation

Weight display mounting bracket included for wall/desk mount.


6800 PORTABLE BARIATRIC FLOOR SCALE

Ideal for any clinical application, this scale is low profile for patient comfort when stepping on and off the scale. With a built-in carrying handle, Detecto's model 6800 scale is easy to transport from office to office. Utilizing Detecto's model 750 weight display and a universal bracket (included), the indicator may be mounted on a wall or desk for easy viewing. Powered by "AA" batteries, this scale is truly portable when power outlets are not accessible. It comes with a 10-foot-long cable from the indicator to the scale.

FEATURES:

- Highly accurate digital weight display with Body Mass Index calculation and easy-to-read LCD.
- Lightweight design, battery power, and carrying handles make it ultra-portable for convenient patient weighing wherever needed.
- Low-profile platform provides easy access.
- Skid-proof and soft, but durable black rubber mat for patient ease and comfort.

COMPACT DESIGN

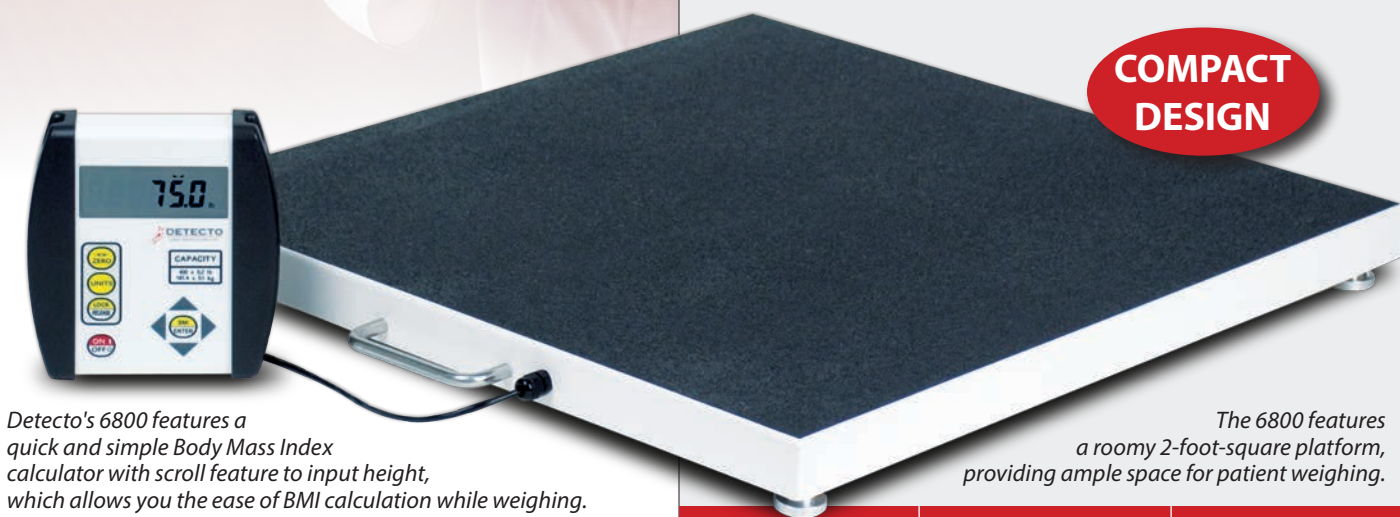

Detecto's 6800 features a quick and simple Body Mass Index calculator with scroll feature to input height, which allows you the ease of BMI calculation while weighing.

The 6800 features a roomy 2-foot-square platform, providing ample space for patient weighing.



VERSATILITY

Detecto's 6800 features an integral carrying handle and a lightweight design, so you may transport this battery-operated scale to any office, room, or home. This lets the weighing come to the patient! With an 800 pound capacity, the 6800 may be used for weighing every patient type from child to bariatric.

MODEL:	6800	6800-AC
DESCRIPTION:	Without AC adapter	With AC adapter
CAPACITY:	1,000 lb x 0.2 lb / 450 kg x 0.1 kg	
PLATFORM DIMENSIONS:	24" W x 2.125" H x 24" D / 61 cm W x 5 cm H x 61 cm D	
DISPLAY DIMENSIONS:	6.125" W x 6.75" H x 1.75" D / 16 cm W x 17 cm H x 4 cm D	
WEIGHT DISPLAY:	0.75" high digits - LCD	
SHIPPING WEIGHT:	45 lb / 20.4 kg	
POWER:	AC adapter (included with model 6800-AC) or 6 "AA" size alkaline batteries (not included)	
CONSTRUCTION:	Baked-on, powder-painted steel	

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.



DETECTO
A Division of Cardinal Scale Manufacturing Co.
102 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153
www.Detecto.com

SOLD BY: