

ADA-COMPLIANT WHEELCHAIR SCALE

7150-AC




The 7150-AC ADA-Compliant Wheelchair Scale is designed for facilities that require enhanced wheelchair access, added safety, and NTEP certified dependable medical accuracy.

ADA-focused design support smooth roll-on weighing while maintaining the highly precision expected in clinical environments.

A 40" platform length combined with 1:8 shallow angle ramps improve wheelchair maneuverability. Integrated side guards help keep wheelchairs properly positioned during weighing. The low-profile design supports patient comfort and efficient workflow in healthcare settings.



DETECTO®

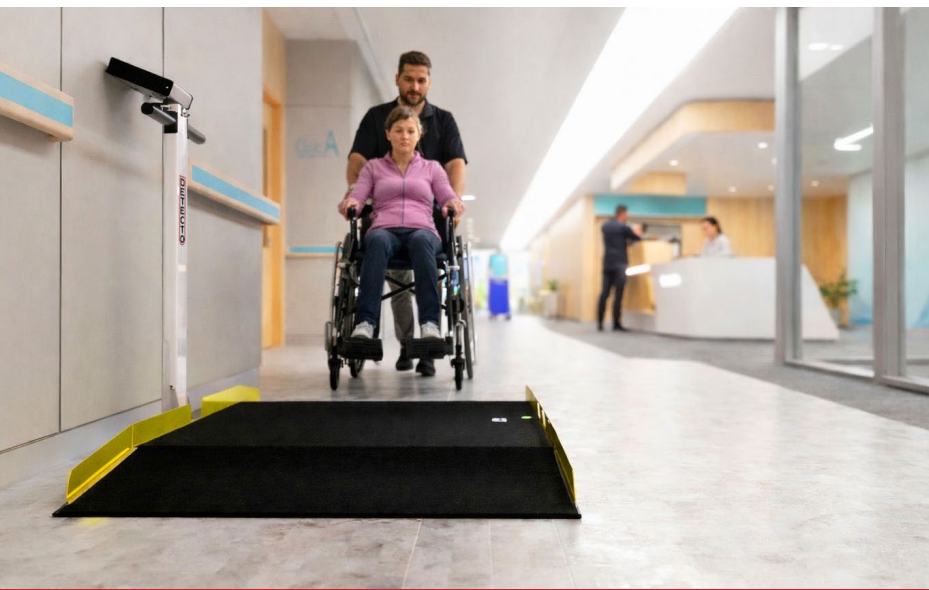

FEATURES

- Side guide rails for additional patient and wheelchair maneuverability
- 1,000 lb / 450 kg Capacity
- 2-way easy-access wheelchair ramps
- Foldable design for easy storage
- Two integral wheels for simple transport
- Easily collapsible column with quick-release lever
- NTEP certified for medical dosage accuracy
- BMI calculation
- EMR/EHR ready
- Low-profile 2" high platform
- OIML EC Class III model available
- Ramps rest on the floor for easy patient use
- Single-pass weighing with push button tare
- Wi-Fi / Bluetooth and Ethernet optional
- Ships 90% assembled

SPECIFICATIONS

- **Capacity:** 1,000 lb x 0.2 lb / 450 kg x 0.1 kg
- **Platform Size:** 40" x 36" x 2.2" (side guards not included in height)
- **Ramp Size:** 36" x 17.4"
- **Folded Depth:** 11.4"
- **Overall Dimensions of Scale:** 75" W x 45" H x 40" D
- **Portability:** Two integral wheels
- **Power:** Model 6800-1045 AC adapter (included) or (6) C batteries (not included)
- **LCD Display Dimensions:** 3.5" W x 3" H
- **LCD Display:** 5-digit, 7-segment, 7/8" high LCD weight digits
- **Keypad:** 19-key, color-coded, sealed membrane
- **Connectivity:** RS232 serial, USB and optional wired ethernet or Wi-Fi
- **Operating Temperature:** 14° to 104°F / -10° to 40°C
- **Net Weight:** 143 lb / 65 kg
- **Shipping Weight:** 190 lb / 86 kg
- **Country of Origin:** USA
- **Resolution:** 5,000 divisions, max.
- **Certifications:** NTEP-certified for dosage accuracy

DETECTO®