

# 8435 SERIES

## DIGITAL BABY SCALES

## DIGITAL SPEED AND ACCURACY



*Easy-to-clean  
polystyrene trays  
cradle babies  
securely on these  
dependable  
pediatric scales.*



DETECTO Models 8435 | 8435KG  
2 models to choose from

### **HEAVY-DUTY SOLID, STEEL SCALE BASE** for patient assurance while weighing.


- Pushbutton tare to remove blanket weights
- Built-in measuring tape: 0 to 21 in x 0.25 in / 0 to 53 cm x 0.1 cm
- Easy-to-clean plastic trays for baby comfort
- Easy-to-read high-contrast LCD
- Automatic weight hold feature for active babies, locks in weight
- Heavy steel base for patient stability
- Durable, long-lasting steel construction
- AC or battery power

# DETECTO DIGITAL BABY SCALES

## DURABLE CONSTRUCTION / CLINICAL QUALITY



The Detecto brand name, known around the world, assures you of a reliable, accurate scale for pediatric weighing. The Detecto models 8435 and 8435KG offer you a choice of pounds or kilograms. Lightweight and attractive, these models feature a heavy-duty, rust-resistant understructure and digital speed and accuracy.

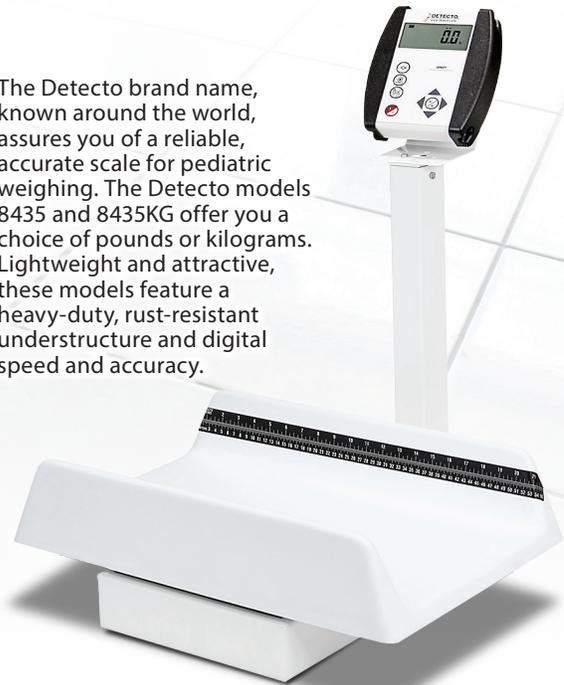

The polystyrene trays are a snap to clean and the scales feature compact footprints. The heavy-duty steel base provides patient assurance while weighing. A built-in measuring tape allows the healthcare professional to measure both weight and height simultaneously.

Detecto's model 8435 and 8435KG feature easy-to-use keypads and extra-large displays perfect for patient weighing and viewing by the attending nurse, physician, or home healthcare professional.

Model	8435	8435KG
Capacity	130 lb x 0.1 lb	59 kg x 0.05 kg
Overall Dimensions	37.5 in H x 22 in W x 21 in D / 95 cm H x 56 cm W x 53 cm D	
Measuring Tape	0 to 21 in x 0.25 in / 0 to 53 cm x 0.1 cm	
Tray Size	22 in L x 14.7 in W x 4.5 in D (56 cm L x 37 cm W x 11 cm D)	
Display	Six-digit, seven segment, 0.7 inch (17.8 mm) high LCD	
Power	Six "AA" Alkaline, Ni-Cad or NiMH batteries (not included) OR an optional 100 to 240V AC 50/60Hz 12V DC 1A wall plug-in UL/CSA listed AC power adapter (Cardinal part number 6800-1045).	
Keyboard	Membrane type with 5 keys and 4 directional arrows	
Connectivity	RS232 serial port	
Construction	Baked-on, powder-painted steel	
Shipping Weight	33 lb (15 kg)	



### CONVENIENT DISPLAY POSITIONING

The weight indicator is conveniently situated over the top of the child, so it's easy to read the weight while the baby is on the scale.



### QUALITY CONSTRUCTION

The 8432-CH has a lock-in weight feature which compensates for typical baby movement. The crisp LCD features 0.75 inch high digits for easy readouts.



### BUILT-IN MEASURING TAPE

The easy-to-read measuring tape allows users to measure both weight and height simultaneously (0 - 21 in/53 cm).

*DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.*



102 E. Daugherty, Webb City, MO 64870 USA  
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153  
[www.Detecto.com](http://www.Detecto.com)

**SOLD BY:**