


BABY SCALES

Digital and mechanical scales for weighing infants to small toddlers

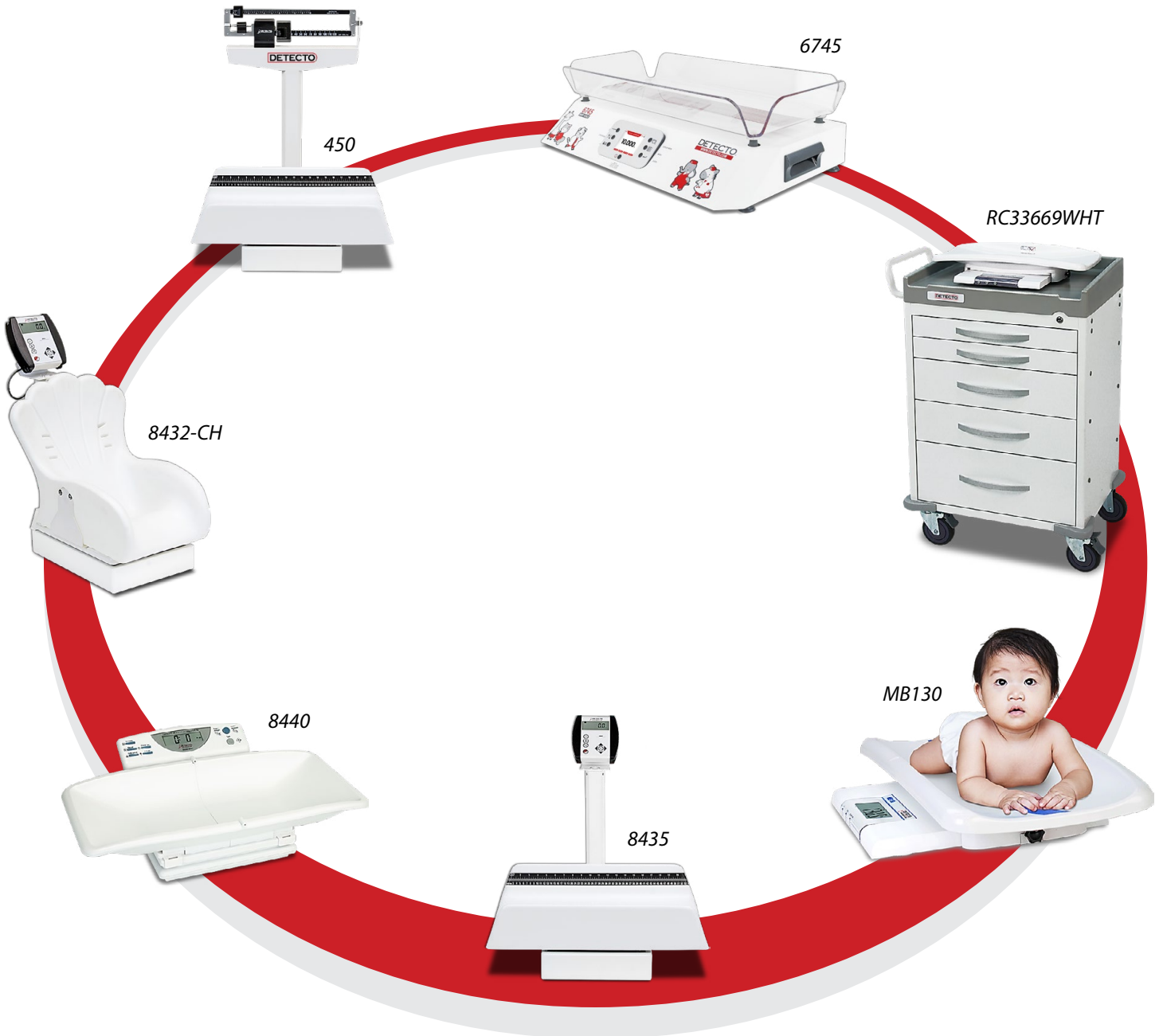

DETECTO[®]

BABY SCALES FOR PEDIATRIC PRACTICES DELIVERY ROOMS NURSERIES AND NICUS

When weighing the most precious of patients, DETECTO provides many different pediatric scales to comfortably weigh and measure babies. DETECTO offers spacious weighing trays to accommodate neonatal and pediatric care units. Weighing trays come with either engraved measuring tapes or stadiometers that can be attached to also gain length measurements. These spacious trays can be horizontal or inclined, while there are economical suspended sling seats and weighing trousers to gain infant weight. There are weigh beam, digital, and suspended configurations depending on pediatrician or clinician needs.

MB130 with DLM

From Neonatal to Pediatric Weighing – Optional Digital Length Measuring Add-on



With dual graduations ranging from neonatal to pediatric care (accurate to 5 grams), DETECTO's MB130 digital scale features a removable tray which automatically converts the scale into a low-profile step-on toddler scale with a maximum capacity of 44 lb / 20 kg. The optional DLM digital length measuring rod can be easily mounted to the scale for an all-digital measuring station. The easy-to-clean plastic tray is 22-inch-wide by 11-inch-deep (56-cm-wide by 29-cm deep). Optional carrying case available.

8440

Standing Toddler or Baby Weighing



Providing accurate and safe weighing of infants and toddlers, DETECTO's 8440 digital scale is ideal for home health and medical clinics. The 8440's weighing tray may be removed for standing toddler weighing and for easier transport. The 1-in / 25-mm high LCD offers pounds and kilograms readouts plus features a built-in clock. The scale's Weight Difference key provides the weight difference between the current weight on the scale and a saved weight in memory. Capacity is 44 lb x 1/2 oz / 20 kg x 10 g.

6745

Premium Quality Baby Scale with .02 oz / 1 g Accuracy



DETECTO's 6745 baby scale was designed for use in the NICU, nursery, delivery room and pediatric practices where high accuracy and dependability is paramount for weighing the smallest of patients. The 6745 features metric lockout, integral measuring tape (up to 22 inch/56 cm), 30 lb x 0.02 oz / 15 kg x 1 g capacity, full-color touchscreen 3.2 in (diagonal) display with 320 x 240 resolution, fun colorful animal graphics, built-in carrying handles, 8-digit numeric patient ID, and roomy 22 in x 9 in / 57 cm x 23 cm weighing tray (flat portion).

8435

Digital Baby Scale with Column-Mounted Display



DETECTO's 8435 digital baby scale features 130 lb x 0.1 lb capacity (metric version also available), pushbutton tare to remove blanket weights, built-in measuring tape (0 to 21 in x 0.25 in / 0 to 53 cm x 0.1 cm), easy-to-clean plastic trays for baby comfort, automatic weight hold feature for active babies, and a heavy steel base for patient stability. The weight indicator is conveniently situated over the top of the child, so it's easy to read the weight while the baby is on the scale.

8432-CH

Inclined Shell-Shaped Seat is Fun for Children



Infants readily accept the inclined seating of DETECTO's 8432-CH series scales. The shell-shaped seat is made of plastic material and fun for children to sit in while weighing. With pound and kilogram readings and heavy-duty steel base and understructure for support and stability, these pediatric scales have a lot to offer. The 8432-CH has a lock-in weight feature which compensates for typical baby movement. The weight indicator is conveniently situated over the top of the seat, so it's easy to read the weight while the baby is seated.

450

Weighbeam Baby Scale Requires No Power



DETECTO's 450 baby scales are available in capacities up to 130 pounds and 65 kilograms with your choice of pounds or kilograms readings. There is no electricity required, making them ideal for remote clinics and pediatric offices. The die-cast weigh beam is positioned high above the scale base to provide for unobstructed viewing of the baby's weight measurement at an easy-to-read height. The baked-on, white powder paint offers a long-lasting dependable finish. The built-in measuring tape measures 0 to 21 in x 0.25 in / 0 to 53 cm x 0.1 cm.

RC33669WHT

Rescue Mobile Baby Scale Cart



The model RC33669WHT baby scale cart perfectly accommodates many of DETECTO's baby scales and digital length measuring rods for a mobile pediatric measurement station throughout your clinical facility. The steel rolling cart with white powder paint features a push handle and two parking locks plus one steering lock for optimal patient safety. The top shelf is made of ABS molded plastic and will securely grip the baby scale's four rubber feet without sliding for patient stability.

Learn More:

<https://detecto.com/product/product-category/Pediatric-Scales>

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

SOLD BY:



102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com