



PORTABLE DIGITAL FLOOR SCALES

DR400C | DR550C



DETECTO DR Series with High-Contrast Backlit Blue LCD

- Compact, low-profile design for space-saving efficiency
- Coiled cable from scale base allows versatile positioning of the indicator
- Slip-resistant, texturized black rubber mat for traction and patient stability
- 550 lb x 0.2 lb / 250 kg x 0.1 kg and 400 lb x 0.5 lb / 180 kg x 0.2 kg models available

PORTABLE

HOME HEALTHCARE SCALES



2.2-in/5.6-cm
high, low-profile
platform for easy
portability

FLAT CONSTRUCTION AND SLIP-PROOF SURFACE FOR COMFORTABLE FOOTHOLD

The portable lightweight platform and large bold backlit blue readouts make this scale ideal for mobile clinics, home healthcare nurses, doctor's offices, general purpose use, and home bathroom weighing. The scale's moveable display lets you place it on a work surface or mounted to the wall, wherever it's most convenient for you to view. An optional black carrying case is available.



DETECTO's model DR400C and DR550C digital portable floor scales will provide you with years of accurate and reliable weighing from the brand you trust in clinical weighing. The scales arrive assembled and ready for use.

- Battery or AC power for your convenience
- Lightweight, mobile weighing platform for easy transport
- Auto Shut-Off after 2 minutes conserves battery power
- Simple four-key operation is easy to use even for novices
- Pushbutton Tare will zero out the load on the platform
- Pounds and kilograms conversion with Units key
- Removeable steel platform that is durable and easy to clean

 **SINCE
1900**

DETECTO is the largest medical scale manufacturer in the United States providing clinical-grade accuracy and dependability for over a century.

Raised buttons and audible beep confirmation of keypresses

The economical DR400C and DR550C are battery powered, so you can take the scales anywhere you need them. The raised buttons are comfortable to press for a pleasant tactile feel, and an audible beep provides a confirmation of your keypresses. The backlit blue LCD with bold one-inch-high digits is viewable in almost any lighting condition.



EXCELLENT FOR USE WHERE SPACE IS LIMITED

Floor/ Tabletop Display

Mobile Positioning

The bracket can temporarily rest on any flat surface for mobile weighing.



Wall Mount Display

Fixed Positioning

Two mounting screws and two plastic drywall anchors are included for mounting the weight display on the wall.



Versatile Display Options

The DR series indicator's versatile drop-in mounting bracket and coiled cable allows the display to be positioned on the floor or tabletop for quick mobile weighing use or you may permanently mount it on the wall using the included screws and plastic drywall anchors.

Battery or AC Power

DETECTO's DR400C and DR550C come with a UL listed 9V 300mA DC adapter (included) or may be powered by 4 AA batteries (not included) for mobile use.

Easy access to 4 AA batteries located conveniently on the back of the indicator.



Four non-skid feet for weighing stability.



PORTABLE

DIGITAL FLOOR SCALES



Includes:

DR series scale with coiled cable to indicator, mounting bracket, DC adapter, operation manual, and two mounting screws and two plastic drywall anchors.

OPTIONAL



Wall-Mount Grab Bar

Stainless steel wall-mount support handles (ordered 1 bar per model)



Model DR-CASE

Carrying Case


Black carrying case with handle and zipper for the models DR400C and DR550C



	DR400C	DR550C
Capacity	400 lb x 0.5 lb / 180 kg x 0.2 kg	550 lb x 0.2 lb / 250 kg x 0.1 kg
Scale Base Platform Dimensions	12.5 in W x 11.9 in D x 2.2 in H / 318 mm W x 303 mm D x 58 mm H	
Weight Indicator Dimensions	6.8 in W x 1.6 in D x 3.2 in H / 172 mm W x 40 mm D x 82 mm H nus bracket)	
Scale Base Platform Construction	Mild steel painted white with textured black rubber mat	
Keys	ON/OFF, HOLD, TARE, lb/kg	
Display Annunciators	Stable weight, Negative weight, Hold, Net, kg, lb, Low battery	
Features	Backlight, low battery indicator, power-saving selectable timing auto shut-off, beeper, printer output	
Display Type	Liquid Crystal Display, 5 digit, 7 segment	
LCD Display Dimensions	3.8 in W x 1.4 in H / 96 mm W x 35 mm H overall LCD size	
LCD Digit Dimensions	1 in H x 0.4 in W / 25 mm H x 10 mm W per digit	
Cable Length	Coiled	34 in / 863 mm
	Fully Extended	90 in / 2,286 mm
Power	9V DC 300mA power supply or (4) AA alkaline batteries (not included)	
Tare	100% of full scale capacity	
Zero	Established on power up routine and maintained by auto-zero circuitry	
Operating Temperature	41° F to 95° F (5° C to 35° C)	
Operating Humidity	35% ~ 90% RH	
Net Weight	10.5 lb / 4.8 kg	
Shipping Weight	12 lb / 5.4 kg	
UPC Code	809161143904	809161143904

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

SOLD BY:

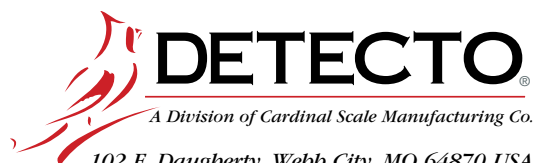

A Division of Cardinal Scale Manufacturing Co.

102 E. Daugherty, Webb City, MO 64870 USA

Pb: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153

www.Detecto.com