

DETECTO®

Free-Standing Portable Height Rod **STADIOMETER**



- Portable, Reliable Height Measurement Accuracy
- 4.5 in – 81 in / 11.5 cm – 205 cm Height Range
- Works in Conjunction with DETECTO's DR400C Scale
- From Pediatric to Adult Measurement
- Height Rod Disassembles Easily into 4 Sections and Stores Conveniently with Head Piece and Base



Optional combo carrying case available: PHR-CASE (holds PHR stadiometer and DR400C floor scale)

Portable Free-Standing **MECHANICAL HEIGHT ROD**

- Measures up to 6'9" / 205 cm
- Two Stabilizers for Maximum Accuracy When Measuring Against a Wall
- Durable and Lightweight All-Plastic Construction
- Measures in Both Inches and Centimeters
- Convenient and Easy to Transport

DETECTO's unique PHR portable height rod is simple and easy to set up. There's no wall fastening necessary and it may be positioned against a wall or completely free standing. Two spacers provide additional stability when placed against a wall. The large floor plate ensures complete stability and permits the height rod to be used anywhere.

Ideal for:

- School Health
- Pediatricians
- General Practitioners
- Mobile Clinics
- Visiting Nurses
- Hospitals


The PHR may be used in conjunction with DETECTO's DR400C digital scale with 12 x 12-inch low-profile platform, 400 lb x 0.5 lb capacity, and one-inch-high LCD readouts. The DR400C scale is powered by 6 AA batteries, so using it together with the PHR height rod provides you with a portable patient measuring station.



DETECTO[®]

Mobile Measurement Convenience AT AN **ECONOMICAL PRICE**

Height measurements are readable on both sides of the PHR height rod.

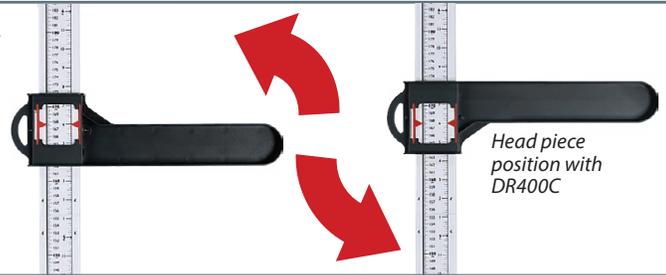

Reversible Head Piece to Accommodate the DR400C Scale

DETECTO's model PHR stadiometer features a unique flipping head piece that lets you add a DETECTO DR400C digital floor scale at the base to measure weight and height simultaneously. If you use the PHR without a DR400C scale at the base, simply flip the head piece upside down to convert it for accurate height readings only. The head piece is sized exactly to compensate for the specific height of the scale platform when used with the DR400C low-profile scale.

Smooth-Gliding Head Piece

Adjusting the smooth-gliding measuring head piece is easy due to the non-warping plastic construction of the height rod. Measurements are readable from either side of the height rod in bold, clearly-visible numbering with bright-red pointer for optimal patient measurement.

PORTABLE CONVENIENCE

DETECTO's optional PHR-CASE carrying bag with zipper and adjustable shoulder strap for convenient transport. Made of black ballistic nylon, the PHR-CASE features pockets fitted exactly for DETECTO's PHR and DR400C.



Optional CARRYING CASE



The carrying case is water and dust resistant providing tailor-made safe storage and compact transport capabilities for these measurement devices. The exclusive PHR-CASE is an indispensable accessory when you purchase DETECTO's PHR and DR400C.

The measuring rod can be taken apart in just a few quick steps and conveniently carried anywhere in the optional carrying case with shoulder strap (model PHR-CASE). The carrying case is cleverly designed with internal pockets that secure the height rod column on one side and the DR400C digital scale and the rest of the stadiometer on the other side.

Mobile Clinical MEASUREMENT



The collapsible height measuring rod of the stadiometer breaks down into four pieces, so in just a few steps it can be attached to the stable platform. When disassembled, the stadiometer is easy to transport.


The parts snap together perfectly for transport ease. Everything fits conveniently into DETECTO's optional PHR-CASE for mobile use. All of these features make this stadiometer the ideal tool for anyone who needs a reliable measuring device for use at different locations.

Model	PHR
Height Range	4.5 in - 81 in / 11.5 cm - 205 cm (visible on both sides of height rod)
Graduation	0.125 in / 0.1 cm
Construction	Plastic
Number of Pieces	8 total
Product Dimensions	13.0 in W x 19.8 in D x 83.5 in H / 33.0 cm W x 50.3 cm D x 212.1 cm H
Product Weight	8.5 lb / 3.8 kg
Storage Dimensions	15 in W x 23.5 in D x 3.8 in H / 38.1 cm W x 59.7 cm D x 9.7 cm H
Shipping Dimensions	14 in W x 23.5 in D x 6.25 in H / 35.5 cm W x 59.7 cm D x 15.9 cm H
Shipping Weight	12 lb / 5 kg
UPC Code	809161199208



Optional
PHR-CASE

Model	PHR-CASE
Construction	Black ballistic nylon with carrying shoulder strap
Storage Content	Can hold both PHR height rod and DR400C scale (pockets included for each)
Product Dimensions	25.0 in W x 6.3 in D x 16.0 in H / 64.3 cm W x 16.0 cm D x 40.6 cm H
Product Weight	4.0 lb / 1.8 kg
Shipping Dimensions	17 in W x 13 in D x 10 in H / 43.1 cm W x 33 cm D x 25.4 cm H
Shipping Weight	8 lb / 3.6 kg
UPC Code	809161199307

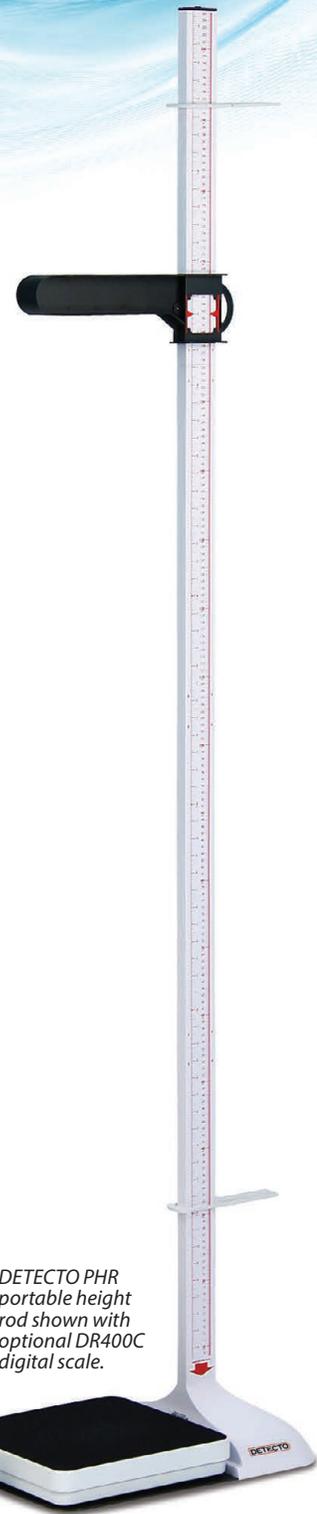

DETECTO PHR portable height rod shown with optional DR400C digital scale.

DETECTO reserves the right to improve, enhance or modify features and specifications without prior notice.

SOLD BY:



102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com