

# DETECTO®

## SOLO-RI LOW-PROFILE SCALE WITH REMOTE INDICATOR



**SOLO**  
DIGITAL CLINICAL SCALES

- Low-profile scale base with 7.25 ft / 2.2 m cable to remote display
- Display may be wall-mounted or set on countertop (bracket included)
- Body Mass Index (BMI) calculation to determine health status
- 14.2 in W x 15.8 in D x 2.4 in H / 36.0 cm W x 40.2 cm D x 6.1 cm H platform size
- 550 lb x 0.2 lb / 250 kg x 0.1 kg capacity

# MEDICAL-GRADE ACCURACY:

## WEIGHT, HEIGHT, AND BODY MASS INDEX

### SOLO-RI Portable Clinical Scales with Space-Saving Compact Design for Flexible Use

- 550 lb x 0.2 lb / 250 kg x 0.1 kg capacity
- Weight, Height, and Body Mass Index
- Low-profile platform (only 2.4 in / 6.1 cm high)
- 7.25 ft (2.2 m) cable from scale base to display to allow placing display wherever it's easiest to read
- Mother/Baby mode factors the weight of infants and toddlers held by an adult
- Wall or desk mount bracket for the display (included)
- Manual input height entry from 2 ft 0.0 in - 7 ft 10.5 in by 0.5 in (60.0 cm - 240.0 cm by 1cm)
- Audible beep confirmation (short beep when a key is pressed, and a long beep for when weight is stable and locked)
- Selectable Auto Shut-Off feature to conserve batteries (up to 9 minutes)
- Auto Weight Lock feature on stable weight
- 6 AA batteries included. SOLO-RI-AC model comes with AC adapter
- Included carpet feet improve weighing accuracy on carpeting

The SOLO-RI provides BMI (Body Mass Index) calculation by conveniently selecting your height while weighing with the easy-to-use scroll arrow keys.



DETECTO's SOLO-RI is portable and light weight for easily changing weighing locations based on patient needs.

**SOLO**  
DIGITAL CLINICAL SCALES



The 7.25 ft/2.2 m cable from the scale base to the indicator allows you to set the display wherever it's most convenient.




The SOLO-RI's remote display may be positioned on countertops, desks, or mounted on the wall for flexible weighing use.

Countertop bracket with rubber feet included for level surfaces



**DETECTO**

*The Name Brand Preferred  
by Clinical Professionals  
Around the World*

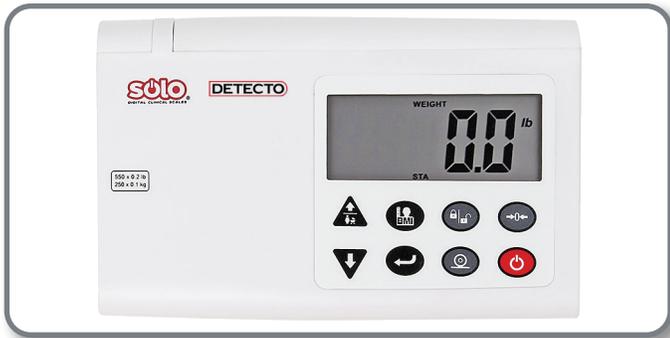

### Stylish Indicator

The SOLO-RI indicator features 8 simple buttons that are easily identifiable in any language for use around the world.



### OneWeigh Units Locking

OneWeigh units locking into Pounds or Kilograms avoids dosage errors based on incorrect units.



### Wall-Mounting Bracket

A multifunctional bracket is included to mount the display on a wall, or the bracket (and included feet) can be used to place it on countertop.



### Weight Hold Key

The SOLO-RI's Hold/Release key lets you hold a patient's weight onscreen until you can record it in their recordkeeping.



### EMR/EHR-Ready

The scale's USB port allows you to output measurement data to EMR/EHR software by pressing the Print key.



### Wide, Flat Platform

The oversized 14 x 15 inch platform with ribbed rubber mat is only 2.4 inches high making it easy for patients to step on/off.



### Everything Included

The SOLO-RI comes ready for weighing right out of the box with included batteries, manual, wall-mount and desktop bracket, screws, and carpet feet.



### Battery Power

No power outlet available, no problem! The SOLO-RI comes with 6 AA batteries for weighing anywhere (model SOLO-RI-AC comes with AC adapter).



Model  
SOLO

DETECTO's  
SOLO-RI is also  
available with an  
eye-level column  
as model SOLO



Model  
SOLO-RI

Optional  
carrying case  
with internal  
Velcro straps  
available (model  
PRODOC-CASE)



Model	SOLO-RI	SOLO-RI-AC
Capacity	550 lb x 0.2 lb / 250 kg x 0.1 kg	
Platform Size	14.2 in W x 15.8 in D x 2.4 in H / 36.0 cm W x 40.2 cm D x 6.1 cm H	
Remote Display Size (with mounting bracket)	7.8 in W x 5.1 in D x 3.0 in H / 19.8 cm W x 12.9 cm D x 7.7 cm H	
Weight Units	Pounds or Kilograms (selectable upon initialization and then locked in)	
Power	(6) AA size Alkaline batteries (included)	(6) AA size Alkaline batteries (included) and medical device 9V AC/DC wall plug-in adapter (DETECTO part number PD-ACKIT, includes UK plug, EU plug, and USA plug)
Display	Five-digit, seven-segment, 1.0-inch (25-mm) high LCD digits	
Operating Temperature	50 to 104 °F / +10 to +40 °C	
Function Keys	ON/OFF, ZERO, BMI HEIGHT, ENTER, LOCK/RELEASE, PRINT, UP ARROW/Mother and Baby, DOWN ARROW	
Annunciators	ZERO, LOCK, STA, HOLD, LB, KG, CM, FT/IN, HEIGHT, WEIGHT, BMI, AM, PM	
Interface	USB (ASCII Protocol)	
Height Units	Inches or centimeters	
Manual Input Height	2 ft 0.0 in – 7 ft 10.5 in by 0.5 in / 60.0 cm - 240.0 cm by 1 cm	
Minimum Weighing	11 lb / 5 kg	
Selectable Auto Shut-Off	Up to 9 minutes	
Construction	base – black ABS; mat – black PVC	
Product Weight	11 lb / 5 kg	12 lb / 5.5 kg
Shipping Weight	14 lb / 6.4 kg	15 lb / 6.9 kg
Package Dimensions	24 in W x 17 in D x 7 in H / 61 cm W x 43 cm D x 18 cm H	
UPC Code	809161330304	809161330700

*Detecto Scale reserves the right to improve, enhance, or  
modify features and specifications without prior notice.*

**SOLD BY:**



102 East Daugherty St., Webb City, MO 64870 USA  
Toll-Free: (800) 641-2008 E-mail: [detecto@cardet.com](mailto:detecto@cardet.com)  
[www.Detecto.com](http://www.Detecto.com)