





## IN-FLOOR DIALYSIS SCALES

DETECTO's Solace series in-floor dialysis scales are flush-mounted platform scales that feature stainless steel decks and include a pit frame and trim ring. The flush surface eliminates tripping hazards and hallway accessibility issues.

The Solace series is available with or without a handrail and all models include a recessed-mount indicator and printer combination. These in-floor scales easily accommodate patients with walkers, canes, wheelchairs, or seated in dialysis chairs.



The ID models are designed to be mounted in a pit flush with the floor. A tubular handrail with padded handgrip is available on some models for walk-on patients.

0.0

## **DIALYSIS + BARIATRIC WEIGHING**

- 1,000 lb x 0.2 lb / 450 kg x 0.1 kg capacity
- 8 standard models to choose from
- Standard platform sizes:
  3 x 3 ft, 4 x 3 ft, 4 x 4 ft,
  and 6 x 4 ft
- Handrails available for patient stability while weighing
- Indicator/printer, pit frame, scale, and trim ring included
- Custom sizes and configurations are available



DETECTO's Solace series in-floor scales are covered with durable stainless steel to provide years of maintenance-free use. The built-in thermal printer makes it easy to print out a weight record for patient recordkeeping.

- Multi-color, touchscreen LCD display
- Built-in thermal ticket printer (included)
- Displays weight, height, and BMI onscreen simultaneously
- Large 0.9-inch-high, clinicalblue LCD weight digits
- Push button or manual entry of tare weights for wheelchairs and stretchers
- Automatically locks on weight temporarily for a few seconds in order to easily view and record measurement
- Menu button provides quick access to user-friendly menu settings for scale customization
- Standard 110V AC power 110V – 240V AC universal switching power supply

ID: 5432 GROSS: 489.0 1b TARE: 100.0 1b NET: 389.0 1b

ID: 8976 GROSS: 218.2 lb TARE: 50.0 lb NET: 168.2 lb

The Solace series in-floor scale's ticket printer can print Gross, Tare, Net and Patient IDs for record keeping.



The Solace series is flush mounted with the floor and includes the pit frame and trim ring, plus an interface cable from the scale to indicator.



DETECTO's Solace series in-floor scales are ideal for dialysis clinics, emergency rooms, and clinical facilities where fast, efficient patient measurement is necessary.





Model

**Capacity** 


**Platform Size** 

Handrail Power

**UPC** Code


809161203707

*In-floor scales comes with model* 855RMP digital weight indicator/ printer combo, pit frame, scale, stainless steel deck, and trim ring. Custom sizes and configurations are available.

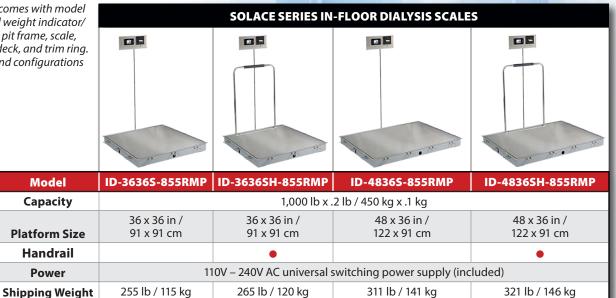

809161203905

809161204001

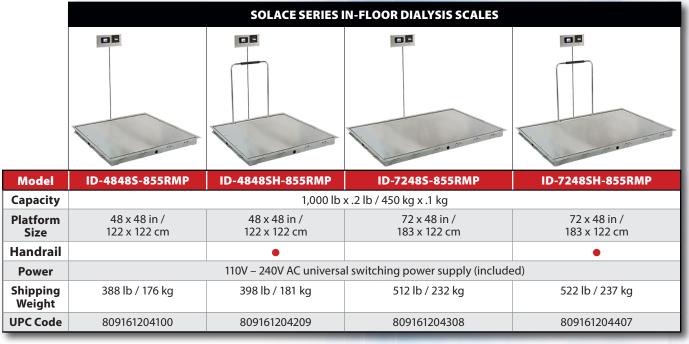

809161203806

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.



203 E. Daugherty, Webb City, MO 64870 USA Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153

Mobile: m.detecto.com

## **SOLD BY:**