

# DETECTO®

## STAINLESS STEEL WEIGHBEAM PHYSICIAN SCALES



Model 439S



Model 337S

**Ideal for Locker Rooms and High-Humidity Areas**

Timeless, attractive appearance constructed in durable stainless steel for optimal weighing accuracy

**STAINLESS  
STEEL** 

# STAINLESS STEEL EYE-LEVEL HEALTHCARE SCALES



**STAINLESS  
STEEL**

Model 439S

## PROTECTION

Rust-resistant stainless steel protects the scale in high-humidity areas like locker rooms.

## HYGIENE

Time-proven benefits of stainless steel: cleans easily and completely for optimal hygiene.

## BEAUTY

Attractive stainless steel scale fits well in all types of environments with timeless appearance.

**DETECTO is the brand name that means Scales in the healthcare industry since 1900.**

Enjoy the classic elegance and added benefits of stainless steel in your DETECTO eye-level weighbeam physician scale. Water damage from mopping around the scale base and high-humidity areas around swimming pools are a thing of the past with DETECTO's stainless steel healthcare scales.

- Trusted by physicians since 1900
- Available in lb, lb/kg, and kg models
- 59 in/148 cm high for eye-level convenience
- Stainless steel construction for long-lasting quality
- Available with or without a height rod
- Die-cast beam can be easily read from either side

Model 439S






Model 3375

# ELEGANT ATTRACTIVE ACCURATE PROTECTION

## MEDICAL WEIGHING PRECISION SINCE 1900

DETECTO's stainless steel eye-level weighbeam scales may be used in gyms, locker rooms, spas, hospitals, clinics, doctor's offices, and even home use.

- 6 different models available in stainless steel
- Capacity of 400 lb/175 kg, 400 lb, or 180 kg
- Simple to operate – no power required
- Easy to read due to eye-level display
- Sturdy construction with stainless steel frame
- Height measurement range: 30-78 in / 76-200 cm
- Height graduations: 0.25 in / 0.5 cm
- Anti-tip column design

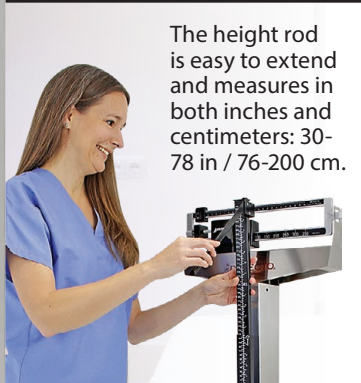

The height rod is easy to extend and measures in both inches and centimeters: 30-78 in / 76-200 cm.



The anti-tip sturdy stainless steel scale base is built rock solid for patient support and stability while weighing and measuring.



The die-cast weighbeam is easy to read from both sides and makes weighing fast and simple (model 3375 shown here).



The slip-resistant black plastic platform cover lifts off easily for cleaning.



Quick and easy calibration via the adjustment screw on the side of the weighbeam.

# SPECIFY "DETECTO"

IT'S THE NAME THAT MEANS SCALES IN THE HEALTHCARE INDUSTRY.



Model	439S	2391S	437S	2371S	339S	337S
Capacity	400 lb x 4 oz	180 kg x 100 g	400 lb x 4 oz	180 kg x 100 g	400 lb x 4 oz / 175 kg x 100 g	400 lb x 4 oz / 175 kg x 100 g
Height Rod	●	●			●	
Shipping Weight	38 lb / 17 kg	38 lb / 17 kg	35 lb / 15 kg	35 lb / 15 kg	38 lb / 17 kg	35 lb / 15 kg
UPC Code	809161132304	809161132403	809161132502	809161132601	809161132700	809161132809

**STAINLESS STEEL**

- Height Measurement Range: 30-78 in / 76-200 cm
- Height Graduations: 0.25 in / 0.5 cm
- Platform Size: 14.75 in x 10.75 in / 37.5 cm x 27.3 cm
- Overall Width: 17.5 in / 44.5 cm
- Overall Height: 58.25 in / 148 cm
- Overall Depth: 20.5 in / 52 cm
- Construction: Stainless steel
- Country of Origin: USA

DETECTO's high-quality weighbeam medical scale is globally recognized as the leader in the mechanical physician scale market.

The innovatively-designed scale features readings from either side of the weighbeam and is imprinted on both sides, allowing the nurse and patient to both view weight simultaneously. Pounds, kilograms, or pounds/kilograms are all available on the easy-to-read beam.

The dual-reading height rod measures in inches and centimeters at a glance (30-78 in / 76-200 cm) with 1/4 inch and 1/2 cm graduations. The slip-resistant black platform cover may be removed for easy cleaning.

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

When the patient steps on the scale, simply move the large poise and then the small poise until the weighbeam balances.

It's that simple!



# DETECTO®

102 East Daugherty St., Webb City, MO 64870 USA  
 Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com  
 www.Detecto.com

**SOLD BY:**