



MEDICAL WEIGHING PRECISION SINCE 1900



DETECTO is the name brand that means Scales in the healthcare industry since 1900. DETECTO's high-quality weighbeam medical scale is globally recognized as the leader in the mechanical physician scale market.

The innovatively-designed scale features readings from either side of the weighbeam and is imprinted on both sides, allowing the nurse and patient to both view weight simultaneously. Pounds, kilograms, or pounds/kilograms are all available on the easy-to-read beam.

The dual-reading height rod measures in inches and centimeters at a glance (24-84 inch / 60-212 cm). The slipresistant black platform cover may be removed for easy cleaning.

DETECTO



When the patient steps on the scale, simply move the large poise and then the small poise until the weighbeam balances. It's that simple!

- 13 different models available
- Capacity of 440 lb/200 kg, 450 lb, or 200 kg
- Simple to operate no power required
- Easy to read due to eye-level display
- Sturdy construction with mild steel frame and powder coat paint finish
- Height measurement range: 24-84 inch / 60-212 cm
- Height graduation: 0.125 inch / 0.1 cm



The slip-resistant steel base and anti-tip column design allows patients to be safely measured and weighed.

IT'S THE NAME THAT MEANS SCALES IN THE HEALTHCARE INDUSTRY.





The die-cast weighbeam is easy to read from both sides and makes weighing fast and simple.

Both the nurse and patient can view the weight at the same time since it is imprinted on both sides.

DETECTO's eye-level weighbeam scales can be used in hospitals, clinics, doctor's offices, gyms, and even home use.



HEIGHT ROD

The height rod is easy to extend and measures in inches and centimeters: 24-84 inch / 60-212 cm



WHEELS

Two rear wheels on some models allow easy transport. The light weight of the DETECTO physician scales combined with wheels makes transport simple.



An important aid for elderly or unsteady patients. The handpost is bolted to the platform, so the patient can steady themselves on it while weighing.




PLATFORM COVER

Slip-resistant platform cover that may be removed for easy cleaning



EYE-LEVEL WEIGHBEAM PHYSICIAN SCALES



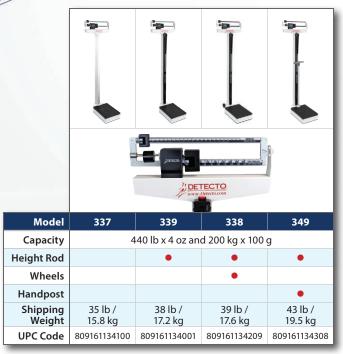

Height Measurement Range: 24-84 in (60 - 212 cm)

Height Graduation: 0.125 in (0.1 cm)

Platform Size: 14.75 in x 10.75 in (37.5 x 27.3 cm)

Overall Width: 17.5 in / 44.5 cm
Overall Height: 58.25 in / 148 cm

Overall Depth (with wheels): 23.25 in / 59 cm

Construction: Mild steel with powder coat paint finish

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

A Division of Cardinal Scale Manufacturing Co.

102 E. Daugherty, Webb City, MO 64870 USA

Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001 www.Detecto.com

SOLD BY: