

# WHEELCHAIR SCALES

Stationary, portable, and wall-mounted scales for wheelchair bound patients that can be used in a variety of healthcare settings



**DETECTO**<sup>®</sup>

# WHEELCHAIR SCALES FOR STANDARD SIZE AND BARIATRIC MOTORIZED WHEELCHAIRS

DETECTO manufactures a complete line of high-quality wheelchair scales to aid medical personnel in weighing patients in wheelchairs. DETECTO offers stationary, portable, and wall-mounted wheelchair scales that can be utilized in a wide variety of healthcare settings. These scales have extra-grip matting for patient safety and dual ramps for easy access from multiple directions. All models are built to use with any standard-size wheelchairs (and some accommodate power wheelchairs) to provide the utmost in comfort, long-lasting durability, and easy portability for wheelchair-bound patients.

## 6550

**Column-Mounted Weight Display with Collapsible Column for Transport**

DETECTO's 6550 portable wheelchair scales offer time-saving and patient comfort benefits. The 6550 features 1,000 lb / 450 kg capacity, Wi-Fi/Bluetooth on some models, vertical fold-up storage, two-way built-in ramps for accessibility from both sides, and 32 in x 36 in / 81 cm x 91 cm platform. With only one hand you can quickly slide the release lever and fold down the 6550's column down for easy transport or vertical storage.



## 6400

**Quick and Easy Mobility for Transport or Vertical Storage – Basic Indicator**


DETECTO's 6400 portable, low-profile wheelchair scale features a high 1,000-lb/450-kg capacity, two-way easy-access ramps, flexible cable remote display for tabletop or wall-mount viewing, 2.2-in/5.6-cm thin platform for easy accessibility, and 32 in x 36 in / 81 cm x 91 cm platform size. Since the display is connected to the scale by a long coiled cable that extends up to 14 ft / 4.2 m, it can be positioned where the healthcare professional can most easily read the results.



## BRW1000

**Larger-Size 32 in x 40 in Platform for Extra-Wide Bariatric Wheelchairs**

DETECTO's BRW1000 is designed to accommodate extra-wide bariatric wheelchairs with its large 32 in x 40 in x 2.2 in / 81 cm x 101 cm x 5.6 cm platform. Since the display is connected to the scale by a long coiled cable that extends up to 14 ft / 4.2 m, it can be positioned where the healthcare professional can most easily read the results. Models are available with Wi-Fi/Bluetooth and also Welch Allyn connectivity. The BRW1000 can accommodate motorized wheelchairs.

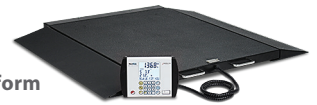

## 6570

**3-in-1 Seated, Standing, and Wheelchair Weighing with Fold-up Seat**

DETECTO's 6570 is a 3-in-1 scale providing seated, standing, and wheelchair weighing. The elegant, curved handrail provides patient stability while weighing standing. The seat folds down for seated weighing and folds back up for wheelchair weighing. By sliding the release levers, the columns easily fold down and lock for fast-and-easy transport and vertical storage. Wi-Fi/Bluetooth is available on some models for wireless EMR/EHR connectivity.



## 7550

**Fold-Up, Wall-Mount Scale Blends Into Any Healthcare Environment**

You can install the 7550 wheelchair scale right in high-traffic areas and fold it up when done weighing. The wall-mount scale with Wi-Fi and Bluetooth stores easily out of the way and aesthetically blends into any environment for a more efficient workflow. The slim 7-inch / 18-cm profile takes up minimal space when installed in a medical facility hallway, clinical office, or patient room. The fold-up non-weighing mode allows you to store the scale in high-traffic areas without using up valuable floor space.



### Learn More:

<https://detecto.com/product/product-category/WheelChair-Scales>

*DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.*

**SOLD BY:**

# DETECTO®

102 East Daugherty St., Webb City, MO 64870 USA  
Toll-Free: (800) 641-2008 E-mail: [detecto@cardet.com](mailto:detecto@cardet.com)  
[www.Detecto.com](http://www.Detecto.com)