

DETECTO®

icon®

DIGITAL CLINICAL SCALES


STATE-OF-THE-ART EYE-LEVEL DIGITAL CLINICAL SCALES

- Extra-high 800-lb capacity
- Sonar touchless height rod for ultimate speed and accuracy
- Ultra-thin platform only 1.5 inches high
- Multi-color touchscreen LCD display
- Displays weight, height, and BMI onscreen simultaneously
- Units locking into pounds or kilograms
- NTEP and OIML certified accuracy for patient dosing



WEIGHT, HEIGHT + BMI INSTANT READOUTS WITHOUT PUSHING A SINGLE BUTTON

The icon's powerful sonar touchless height rod combined with an all-in-one color display means the patient simply steps onto the scale platform and within a second your full measurements are all displayed onscreen instantly without having to press a single key for the ultimate in speed, hygiene, and accuracy.

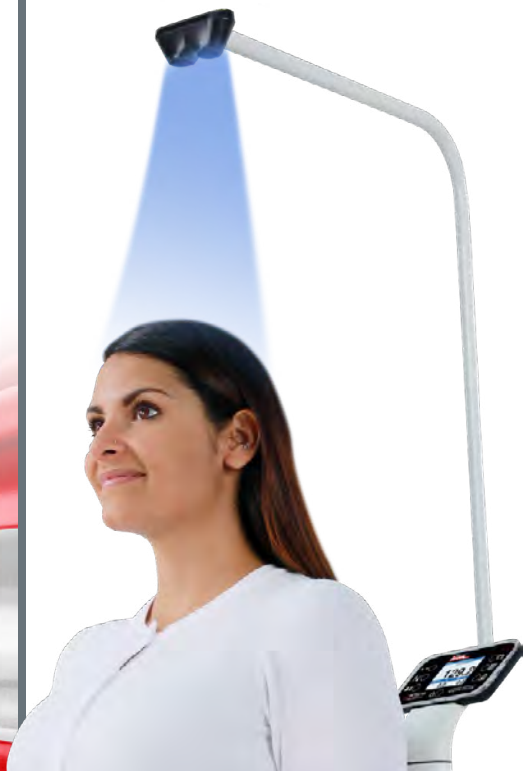



MINIMAL STEP

The icon's virtually-flat base is only 1.5 inches high, offering the ultimate in patient ease of use and comfort while weighing. Plus, the extra-wide 17 inch wide by 17 inch deep platform size allows bariatric weighing.

THE NEXT EVOLUTION OF CLINICAL PATIENT MEASUREMENT HAS ARRIVED!

STATE-OF-THE-ART MEASUREMENT TECHNOLOGY THAT IS ICONIC

Every so often in history, the evolution of product technology takes a monumental step forward. DETECTO has been pioneering physical measurement since 1900 and along the way studying intensively what comprises an incredible healthcare weighing experience. As result of over a century of research, we introduce you to clinical weighing's new icon®.

DETECTO's icon® eye-level digital clinical scale with sonar height rod has an extra high capacity of 800 lb x 0.1 lb / 400 kg x 0.05 kg. The low-profile 17 in W x 17 in D platform is nearly flat at only 1.5 inches high. The icon® includes a non-medical grade standard AC adapter or can be powered by 12 AA batteries (not included). The touchless sonar height rod allows you to simply step on the scale and within seconds weight, height, and Body Mass Index will all be displayed onscreen without touching a single key. The Mother/Baby function factors the weight of infants and toddlers held by an adult.



DETECTO's icon® scales measure weight and height efficiently on patients of all ages and sizes.



CONNECTIVITY PORTS STANDARD IN THE BASE

DETECTO's icon® represents the next generation of digital clinical measurement. There is nothing equivalent to the icon® on the market!

- Sleek, stylish product design
- Eye-level weight display convenience
- Sturdy column construction
- 3 different standard icon® models available to choose from
- Mother/Baby function
- AC adapter included



MULTI-COLOR TOUCHSCREEN DISPLAY

The icon® features a high-tech, smart-phone-style indicator for easy utilization with 7 simple buttons and a multi-color touchscreen LCD display for large readouts.

- 0.9-inch-high, clinical-blue LCD weight digits
- Multi-color touchscreen LCD display
- Displays weight, height, and Body Mass Index onscreen simultaneously
- Automatically locks on weight temporarily for a few seconds in order to easily view and record measurement
- Menu button provides quick access to user-friendly menu settings for scale customization
- Battery power level indication

VERSATILE CONNECTIVITY:

- 1 USB-B port
- HL7 IEEE 11073 compliant (standard)
- Wi-Fi/Bluetooth for wireless EMR/EHR
- Welch Allyn® device connectivity available: CVSM, CSM, and LXI (see models chart on back page)

PATIENT IDs



The icon's brilliant, multi-color touchscreen display and intuitive user-prompted navigation menus provide easy entry of items such as 8-digit numeric patient IDs for EMR/EHR.

AC ADAPTER



DETECTO's icon® clinical scale includes 1 USB-B port, integral bubble level, AC adapter is included standard, Wi-Fi/Bluetooth, and Welch Allyn EMR/EHR connectivity available.

AUTO LOCKING



When a patient steps on the scale, weight, height, and BMI are all shown onscreen and instantly held on the display for 7 seconds automatically in order for the medical professional to view and record the data.

BATTERIES



The platform cover may be removed for easy cleaning and also to provide access to the battery compartment. The icon® may be powered by 12 AA batteries (not included).

STATE-OF-THE-ART EYE-LEVEL DIGITAL CLINICAL SCALES

icon® EYE-LEVEL DIGITAL SCALES WITH SONAR HEIGHT RODS			
Model	ICON	ICON-UWA	ICON-LXI
Capacity	800 lb x 0.1 lb / 400 kg x 0.05 kg		
Height Rod	Touchless sonar height rod		
Height Rod Measures	24" (2') - 86" (7'2") x 0.2 in / 60 cm - 218 cm x 0.5 cm		
DIMENSIONS	Overall Dimensions	17 in W x 18.5 in D x 89.5 in H / 43 cm W x 47 cm D x 227.3 cm H	
	Platform Size	17 in W x 17 in D x 1.5 in H / 43.2 cm W x 43.2 cm D x 3.8 cm H	
	LCD	.9 in / 2.3 cm H LCD digits	
POWER	Batteries	12 AA batteries (not included)	
	AC Adapter (included)	•	•
	Ports	1 RS232 serial port and 1 USB-B port	
	Connectivity	Wi-Fi/Bluetooth	Welch Allyn CVSM/CSM Wi-Fi/Bluetooth
	Net Weight	22 lb / 10 kg	
	Shipping Weight	31 lb / 14 kg	
	Shipping Dimensions	58 in x 20 in x 6 in / 147 cm x 51 cm x 15 cm	
	UPC Code	809161198805	809161199109

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.



Bluetooth™

Wi-Fi/Bluetooth
connectivity included



DETECTO®

102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com

SOLD BY: